

Infitek

MULTIPARAMETER METERS



Infitek



Infitek



Infitek



Infitek



Infitek


Infitek

Infitek Inc.

Address: 522W Riverside Ave Ste N, Spokane, WA 99201, USA
Email: Info@infitek.com

Infitek Co., Ltd.

TEL: +86-531-88982330
FAX: +86-531-88983691
Website: infitek.com
Email: info@infitek.com
Service: service@infitek.com
Address: Room 201A, Building 9 (Block B), No. 2679,
Hechuan Road, Minhang District, Shanghai,P.C.201100

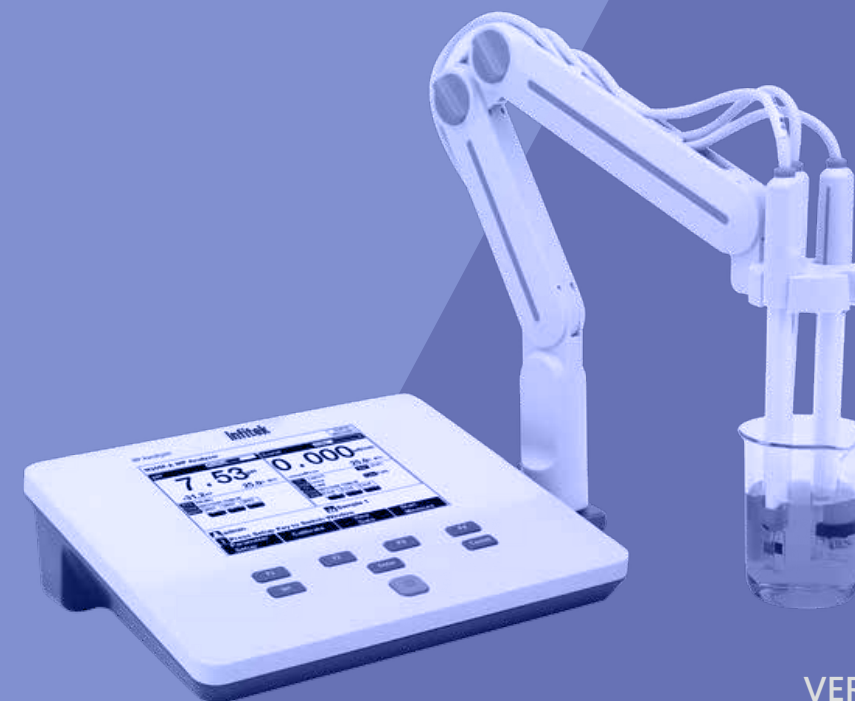

VERSION.26

Benchtop Multiparameter Meter

BEP-M300A

Features

-  High resolution LCD display.
-  Data Storage 500 sets.
-  IP54 waterproof.
-  Support for USB and RS-232 communication.



Description

- The BEP-M300A Benchtop Multi-parameter Analyzer adopts a multi-purpose measurement method, which is suitable for the fields of health disease control, environmental protection, biological agriculture and forestry, water quality analysis, scientific research and institutions of higher learning.

Features

General Features

- High resolution LCD display, 5.7 inches.
- Multi-reading feature allows auto-read, timed-read and continuous-read.
- Automatic/Manual temperature compensation ensures accurate results.
- Auto-hold feature senses and locks the measurement endpoint.
- Data Storage 500 sets (GLP-compliant).
- Support for USB and RS-232 communication.
- Reset feature automatically resumes all settings back to factory default options.
- IP54 waterproof.

pH

- 1-5 points calibration with Standard Recognition.
- Selectable pH buffer groups, including NIST, DIN, GB.
- Automatic electrode diagnosis with pH slope and offset display.

Ion

- 1-5 points calibration.
- Selectable measurement unit, including $\mu\text{g/L}$, mg/L , g/L , mmol/L , PX, etc.
- Multi-measurement modes are supported, including Direct Reading mode, Standard Addition mode, Sample Addition mode and GRAN mode.
- Over 10 methods are built-in, including F^- , Cl^- , Br^- , I^- , NO_3^- , BF_4^- , NH_4^+ , K^+ , Na^+ , Ca^{2+} , Cu^{2+} , Pb^{2+} , Ag^+ and etc., user-defined method is supported.

Conductivity

- 1-3 points calibration with Standard Recognition.
- Settable parameters, including cell constant, temperature compensation coefficient and TDS factor.
- Temperature compensation type (none, linear, pure water).



BEP-M300A Multi-parameter Analyzer 1 Set

- E-301-QC 3 in 1 pH composite electrode 1 Piece
- DJS-IVTC Conductivity electrode 1 Piece
- Standard pH Buffer Solution (50mL) (pH4.01, 7.00, 10.01) 1 Bottle each
- Conductivity standard $1413\mu\text{S/cm}$ (50mL) 1 Bottle
- REX-5 electrode stand 1 Set
- European plug model(QD24-0902000-A1)Power adapter DC 9V (Positive inside and negative outside) 1 Piece
- Universal USB cable 1 Piece
- Dust cover (Key function version) 1 Piece
- Operation quick guide 1 Piece
- Instruction manual 1 Piece

Specifications

Model		BEP-M300A
Parameters		pH/EC/ISE/Temp. (mV/ORP/pX/Resistivity/TDS/Sal.)
pH	Range	-2.00 to 20.00pH
	Resolution	0.1, 0.01 pH
	Accuracy	±0.01 pH
	Calibration Points	Up to 5
	Standard Customization	Yes
	Standard Recognition	NIST, GB and DIN buffers
	Slope Limit	Yes
mV	Range	-2000.0 to 2000.0 mV
	Resolution	0.1 mV
	Accuracy	±0.3 mV or ±0.1%
pX	Range	-2.00 to 20.00
	Resolution	0.1, 0.01 pX
	Accuracy	±0.01 pX
	Calibration Points	Up to 5
ISE	Range	1E-9 to 9.999E9
	Unit	mol/L, mmol/L, g/L, mg/L, µg/L, ppm
	Resolution	Up to 4 significant digits
	Accuracy	±0.5%
	Calibration Points	Up to 5
Conduc- tivity	Range	0.000 µS/cm to 1000 mS/cm
	Resolution	0.001 µS/cm minimum, various with range selection
	Accuracy	±1.0% FS
	Reference Temperature	20, 25 °C
	Calibration Points	Up to 3
	Standard Recognition	84µS/cm, 1413µS/cm, 12.88mS/cm
Resistivity	Range	5.00 Ω·cm~20.00 MΩ·cm
	Resolution	0.01 Ω·cm minimum
	Accuracy	±1.0% FS
TDS	Range	0.00ppm~300ppt
	Resolution	0.01mg/L minimum, various with range selection
	Accuracy	±1.0% FS
Salinity	Range	0~80ppt
	Resolution	0.1ppt
	Accuracy	±2ppt

Model		BEP-M300A
Temperature	Range	-5 to 110 °C, 23 to 230 °F
	Unit	°C, °F
	Resolution	0.1
	Accuracy	±0.2
Measurement	Reading Mode	Auto Read (Fast, Medium, Slow), Timed, Continuous
	Reading Prompts	Reading, Stable, Locked
	Temp. Compensation	ATC, MTC
Data manage- ment	Data Storage	500 results each
	GLP Features	Yes
Inputs	pH Electrode	BNC(Q9)
	DO with Temp. Probe	4-pin aviation connector
	Conductivity with Temp. Probe	5-pin aviation connector
Outputs	USB	USB 2.0 flash memory device, PC
	RS 232	printer
Display options	Backlight	Yes
	Auto Shutdown	1~60 min, off
	IP Rating	IP54
	Date and Time	Yes
General	Electricity	AC Adapter, 100~240V AC input, DC9V output
	Dimensions	242×195× 68 mm
	Weight	900g



Benchtop Multiparameter Meter

BEP-M300F



pH

- ◆ 1-5 points calibration with Standard Recognition.
- ◆ Selectable pH buffer groups, including NIST, DIN, GB.
- ◆ Automatic electrode diagnosis with pH slope and offset display.

Conductivity

- ◆ 1-3 points calibration with Standard Recognition.
- ◆ Settable parameters, including cell constant, temperature compensation coefficient and TDS factor.
- ◆ Temperature compensation type (none, linear, pure water).


Ion

- ◆ 1-5 points calibration.
- ◆ Selectable measurement unit, including $\mu\text{g/L}$, mg/L , g/L , mmol/L , PX , etc.
- ◆ Multi-measurement modes are supported, including Direct Reading mode, Standard Addition mode, Sample Addition mode and GRAN mode.
- ◆ Over 10 methods are built-in, including F^- , Cl^- , Br^- , I^- , NO_3^- , BF_4^- , NH_4^+ , K^+ , Na^+ , Ca^{2+} , Cu^{2+} , Pb^{2+} , Ag^+ and etc., user-defined method is supported.

DO

- ◆ Support for air-saturated water or zero oxygen calibration.
- ◆ Selectable pressure unit, including kPa, mbar, Torr, Atm.
- ◆ Auto barometric pressure compensation
- ◆ Manual Salinity Factor Correction

Features

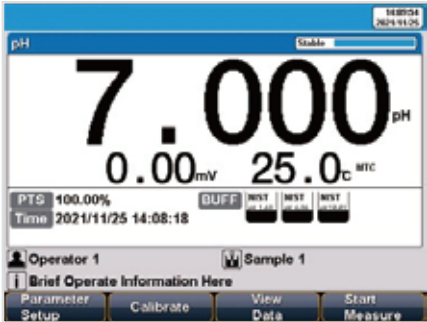

High resolution LCD display, 5.7 inches.
Smart UI design Best user experience, highly efficient measurement



Support for USB and RS-232 communication.
Multiple external devices available PC, printer



Automatic identification standard solution
Fast calibration and more accurate measurement

Description

The BEP-M300F Benchtop Multi-parameter Analyzer adopts a multi-purpose measurement method, which is suitable for the fields of health disease control, environmental protection, biological agriculture and forestry, water quality analysis, scientific research and institutions of higher learning.

General Features

IP54 waterproof.

Multi-reading feature allows auto-read, timed-read and continuous-read.

Automatic/Manual temperature compensation ensures accurate results.

Data Storage 500 sets (GLP-compliant).

Reset feature automatically resumes all settings back to factory default options.

Auto-hold feature senses and locks the measurement endpoint.

Specifications

Model		BEP-M300F
Parameters		pH/EC/ISE/DO/Temp. (mV/ORP/pX/Resistivity/TDS/Sal./DO Saturation)
pH	Range	-2.00 to 20.00pH
	Resolution	0.1, 0.01 pH
	Accuracy	±0.01 pH
	Calibration Points	Up to 5
	Standard Customization	Yes
	Standard Recognition	NIST, GB and DIN buffers
	Slope Limit	Yes
mV	Range	-2000.0 to 2000.0 mV
	Resolution	0.1 mV
	Accuracy	±0.3 mV or ±0.1%
pX	Range	-2.00 to 20.00
	Resolution	0.1, 0.01 pX
	Accuracy	±0.01 pX
	Calibration Points	Up to 5
ISE	Range	1E-9 to 9.999E9
	Unit	mol/L, mmol/L, g/L, mg/L, µg/L, ppm
	Resolution	Up to 4 significant digits
	Accuracy	±0.5%
	Calibration Points	Up to 5
Conductivity	Range	0.000 µS/cm to 1000 mS/cm
	Resolution	0.001 µS/cm minimum,various with range selection
	Accuracy	±1.0% FS
	Reference Temperature	20, 25 °C
	Calibration Points	Up to 3
	Standard Recognition	84µS/cm, 1413µS/cm, 12.88mS/cm
Resistivity	Range	5.00 Ω·cm~20.00 MΩ·cm
	Resolution	0.01 Ω·cm minimum
	Accuracy	±1.0% FS
TDS	Range	0.00ppm~300ppt
	Resolution	0.01mg/L minimum, various with range selection
	Accuracy	±1.0% FS
Salinity	Range	0.0~80.0ppt
	Resolution	0.1ppt
	Accuracy	±2ppt

Model		BEP-M300F
Dissolved Oxygen Concentration	Range	0.00 to 20.00 ppm
	Resolution	0.01ppm
	Accuracy	±0.10ppm
	Calibration Points	Air-saturated water or zero point
	Barometric Compensation	Yes
Saturation	Manual Salinity Factor Correction	Yes
	Range	(0.0 to 200.0)%
	Resolution	0.1%
	Accuracy	±2.0%
Temperature	Range	-5 to 110 °C, 23 to 230 °F
	Unit	°C, °F
	Resolution	0.1
	Accuracy	±0.2
Measurement	Reading Mode	Auto Read (Fast, Medium, Slow), Timed, Continuous
	Reading Prompts	Reading, Stable, Locked
	Temp. Compensation	ATC, MTC
Data manage-ment	Data Storage	500 results each
	GLP Features	Yes
Inputs	pH Electrode	BNC(Q9)
	DO with Temp. Probe	4-pin aviation connector
	Conductivity with Temp. Probe	5-pin aviation connector
Outputs	USB	USB 2.0 flash memory device, PC
	RS 232	printer
Display options	Backlight	Yes
	Auto Shutdown	1~60 min, off
	IP Rating	IP54
	Date and Time	Yes
General	Electricity	AC Adapter, 100~240V AC input, DC9V output
	Dimensions	242×195× 68 mm
	Weight	900g



- BEP-M300F Multi-parameter Analyzer 1 Set**
- E-301-QC 3 in 1 pH composite electrode 1 Piece
 - DJS-IVTC Conductivity electrode 1 Piece
 - DO-958-Q Polarographic DO electrode 1 Piece
 - Standard pH Buffer Solution (50mL) (pH4.01, 7.00, 10.01) 1 Bottle each
 - Conductivity standard 1413µS/cm(50mL) 1 Bottle
 - REX-5 electrode stand 1 Set
 - Power adapter DC 9V(European plug model, QD24-0902000-A1) 1 Piece
 - Universal USB cable 1 Piece
 - Dust cover (Key function version) 1 Piece
 - Operation quick guide 1 Piece
 - Instruction manual 1 Piece

Benchtop Multiparameter Meter

BEP-M500A



pH

- ◆ 1-6 points calibration with Standard Recognition.
- ◆ Selectable pH buffer groups, including DIN, NIST, USA, MERK, JIS, GB.
- ◆ Support for pH standard group customization.
- ◆ Automatic electrode diagnosis with pH slope and offset display.

Conductivity

- ◆ 1-5 points calibration with Standard Recognition.
- ◆ Settable parameters, including cell constant, reference temperature (5/10/15/18/20/25°C), temperature compensation coefficient and TDS factor.
- ◆ Temperature compensation type (none, linear, pure water).

Ion

- ◆ 1-6 points calibration
- ◆ Selectable measurement unit, including $\mu\text{g/L}$, mg/L , g/L , mmol/L , pX , ppm , etc.
- ◆ Multi-measurement modes are supported, including Direct Reading mode, Standard Addition mode, Sample Addition mode and GRAN mode.
- ◆ Over 10 methods are built-in, including F^- , Cl^- , Br^- , I^- , NO_3^- , BF_4^- , NH_4^+ , K^+ , Na^+ , Ca^{2+} , Cu^{2+} , Pb^{2+} , Ag^+ and etc., user-defined method is supported.

Features



Color High contrast LCD touch-screen, 7.0 inches.



Support for USB and RS-232 communication. Multiple external devices available
PC, printer



Automatic identification standard solution
Fast calibration and more accurate measurement

Description

The BEP-M500A Benchtop Multi-parameter Analyzer adopts a multi-purpose measurement method, which is suitable for the fields of health disease control, environmental protection, biological agriculture and forestry, water quality analysis, scientific research and institutions of higher learning.

General Features



IP54 waterproof.



Intelligent operation system provides features including user management, calibration management, method management, data management and etc.



Automatic/Manual temperature compensation ensures accurate results.



Multi-reading feature allows auto-read, timed-read and continuous-read.



Data Storage 1000 sets for each parameter (GLP-compliant).

Data analysis feature helps user review, compare and recalculate results.



Reset feature automatically resumes all settings back to factory default options.

Specifications

Model		BEP-M500A
Parameters		pH/EC/ISE/Temp. (mV/ORP/pX/Resistivity/TDS/Sal.)
pH	Range	-2.000 to 20.000 pH
	Resolution	0.1, 0.01, 0.001pH
	Accuracy	±0.002 pH
	Calibration Points	Up to 6
	Standard Customization	Yes
	Calibration Reminder	Yes
	Standard Recognition	NIST, GB, DIN, USA, MERK, JIS buffers
	Slope Limit	Yes
mV/ORP	Range	-2000.00 to 2000.00 mV
	Resolution	0.1, 0.01 mV
	Accuracy	±0.1 mV or ±0.03% of reading whichever is greater
	EH ORP Mode	Yes
	ORP Calibration Points	1 custom point (relative mV)
pX	Range	-2.000 to 20.000 pX
	Resolution	0.1, 0.01, 0.001 pX
	Accuracy	±0.002 pX
	Calibration Points	Up to 6
ISE	Range	1E-9 to 9.999E9
	Units	mol/L, mmol/L, g/L, mg/L, µg/L, ppm
	Resolution	Up to 4 significant digits
	Accuracy	±0.3%
	Calibration Points	Up to 6
Conductivity	Range	0.000 µS/cm to 2000 mS/cm
	Resolution	0.001 µS/cm minimum,various with range selection
	Accuracy	±0.5% FS
	Reference Temperature	5, 10, 15, 18, 20, 25 °C
	Calibration Points	Up to 5
	Calibration Reminder	Yes
	Standard Recognition	10µS/cm, 84µS/cm, 500µS/cm, 1413µS/cm, 12.88mS/cm
Resistivity	Range	5.00 Ω·cm~100.00 MΩ·cm
	Resolution	0.01 Ω·cm minimum
	Accuracy	±0.5% FS
TDS	Range	0.000~1000ppt
	Resolution	0.001ppm minimum, changed with range
	Accuracy	±0.5%

Model		BEP-M500A
Salinity	Type	NaCl
	Range	0.0~80.0ppt
	Resolution	0.1ppt
	Accuracy	±1ppt
Temperature	Range	-10 to 135 °C, 14 to 275 °F
	Unit	°C, °F
	Resolution	0.1
	Accuracy	±0.1
Measurement	Reading Mode	Auto Read (Fast, Medium, Slow), Timed, Continuous
	Reading Prompts	Reading, Stable, Locked
	Temp. Compensation	ATC, MTC
Data manage- ment	Data Storage	1000 Groups
	GLP Features	Yes
	Log Management	Yes
Inputs	pH Electrode	BNC(Q9)
	DO with Temp. Probe	4-pin aviation connector
	Conductivity with Temp. Probe	5-pin aviation connector
Outputs	USB	USB 2.0 flash memory device, PC
	RS 232	printer, scanner, auto-sampler
Display options	Backlight	Yes
	Auto Shutdown	1~60 min, off
	IP Rating	IP54
	Date and Time	Yes
General	Electricity	AC Adapter, 100~240V AC input, DC24V output
	Dimensions	220×195×68 mm
	Weight	950g



- BEP-M500A Multi-parameter Analyzer 1 Set**
- E-301-QC 3 in 1 pH composite electrode 1 Piece
 - DJS-IVTC Conductivity electrode 1 Piece
 - Standard pH Buffer Solution (50mL) (pH4.01, 7.00, 10.01) 1 Bottle
 - Conductivity standard 1413µS/cm(50mL) 1 Bottle
 - REX-6 electrode stand 1 Set
 - GC65-240270-F Power adapter DC24V 1 Piece
 - Dust cover (Touch screen version) 1 Piece
 - European standard power cord(3VTJ2) 1 Piece
 - Universal USB cable 1 Piece
 - Touch pen (Capacitive screen) 1 Piece
 - Operation quick guide 1 Piece
 - Instruction manual 1 Piece

Benchtop Multiparameter Meter

BEP-M500T



Features

-  Color High contrast LCD touchscreen.
-  Data Storage 1000 sets for each parameter.
-  IP54 waterproof.
-  Support for USB and RS-232 communication.



Description

- The BEP-M500T Benchtop Multi-parameter Analyzer a multi-purpose measurement method, which is suitable for the fields of health disease control, environmental protection, biological agriculture and forestry, water quality analysis, scientific research and institutions of higher learning.

Features

General Features

- Color High contrast LCD touchscreen, 7.0 inches.
- Intelligent operation system provides features including user management, calibration management, method management, data management and etc.
- Multi-reading feature allows auto-read, timed-read and continuous-read.
- Automatic/Manual temperature compensation ensures accurate results.
- Data Storage 1000 sets for each parameter (GLP-compliant).
- Data analysis feature helps user review, compare and recalculate results.
- Support for USB and RS-232 communication.
- Reset feature automatically resumes all settings back to factory default options.
- IP54 waterproof.

pH

- 1-6 points calibration with Standard Recognition.
- Selectable pH buffer groups, including DIN, NIST, USA, MERK, JIS, GB.
- Support for pH standard group customization.
- Automatic electrode diagnosis with pH slope and offset display.

Ion

- 1-6 points calibration
- Selectable measurement unit, including $\mu\text{g/L}$, mg/L , g/L , mmol/L , pX , ppm , etc.
- Multi-measurement modes are supported, including Direct Reading mode, Standard Addition mode, Sample Addition mode and GRAN mode.
- Over 10 methods are built-in, including F^- , Cl^- , Br^- , I^- , NO_3^- , BF_4^- , NH_4^+ , K^+ , Na^+ , Ca^{2+} , Cu^{2+} , Pb^{2+} , Ag^+ and etc., user-defined method is supported.

Conductivity

- 1-5 points calibration with Standard Recognition.
- Settable parameters, including cell constant, reference temperature (5/10/15/18/20/25°C), temperature compensation coefficient and TDS factor.
- Temperature compensation type (none, linear, pure water).

DO

- Support for air-saturated water or zero oxygen calibration.
- Auto barometric pressure compensation.
- Manual Salinity Factor Correction.
- Selectable pressure unit, including kPa , mbar , Torr , Atm .

Specifications

Model			BEP-M500T
Parameters			pH/EC/ISE/DO/Temp. (mV/ORP/pX/Resistivity/TDS/Sal./DO Saturation)
pH	Range	-2.000 to 20.000 pH	
	Resolution	0.1, 0.01, 0.001pH	
	Accuracy	±0.002 pH	
	Calibration Points	Up to 6	
	Standard Customization	Yes	
	Calibration Reminder	Yes	
	Standard Recognition	NIST, GB, DIN, USA, MERK, JIS buffers	
	Slope Limit	Yes	
mV/ORP	Range	-2000.00 to 2000.00 mV	
	Resolution	0.1, 0.01 mV	
	Accuracy	±0.1 mV or ±0.03% ofreading whichever is greater	
	EH ORP Mode	Yes	
	ORP Calibration Points	1 custom point (relative mV)	
pX	Range	-2.000 to 20.000 pX	
	Resolution	0.1, 0.01, 0.001 pX	
	Accuracy	±0.002 pX	
	Calibration Points	Up to 6	
ISE	Range	1E-9 to 9.999E9	
	Units	mol/L, mmol/L, g/L, mg/L, µg/L, ppm	
	Resolution	Up to 4 significant digits	
	Accuracy	±0.3%	
	Calibration Points	Up to 6	
Conduc-tivity	Range	0.000 µS/cm to 2000 mS/cm	
	Resolution	0.001 µS/cm minimum,various with range selection	
	Accuracy	±0.5% FS	
	Reference Temperature	5, 10, 15, 18, 20, 25 °C	
	Calibration Points	Up to 5	
	Calibration Reminder	Yes	
	Standard Recognition	10 µS/cm, 84 µS/cm, 500 µS/cm, 1413 µS/cm, 12.88 mS/cm	
Resistivity	Range	5.00 Ω·cm~100.00 MΩ·cm	
	Resolution	0.01 Ω·cm minimum	
	Accuracy	±0.5% FS	
TDS	Range	0.000~1000 ppt	
	Resolution	0.001 ppm minimum, changed with range	
	Accuracy	±0.5%FS	
Salinity	Type	NaCl	
	Range	0.0~80.0ppt	
	Resolution	0.1ppt	
	Accuracy	±1ppt	

Model			BEP-M500T
Dissolved	Range	0.00 to 99.99 ppm	
Oxygen	Resolution	0.01 ppm	
Concentra-tion -	Accuracy	±0.10ppm, 0.00~50.00ppm ±0.50ppm, 50.0~99.99ppm	
	Calibration Reminder	Yes	
Polaro-graphic	Calibration Points	Air-saturated water or zero point	
	Barometric Compensation	Yes	
	Manual Salinity Factor Correction	Yes	
% Saturation	Range	(0.0 to 600)%	
-Polaro-graphic	Resolution	0.10%	
	Accuracy	±2.0%	
Tempera-ture	Range	-10 to 135 °C, 14 to 275 °F	
	Unit	°C, °F	
	Resolution	0.1	
Measure-ment	Accuracy	±0.1	
	Reading Mode	Auto Read (Fast, Medium, Slow), Timed, Continuous	
	Reading Prompts	Reading, Stable, Locked	
Data manage-ment	Temp. Compensation	ATC, MTC	
	Data Storage	1000 Groups	
	GLP Features	Yes	
Inputs	Log Management	Yes	
	pH Electrode	BNC(Q9)	
	DO with Temp. Probe	4-pin aviation connector	
Outputs	Conductivity with Temp. Probe	5-pin aviation connector	
	USB	USB 2.0 flash memory device, PC	
	RS 232	printer, scanner, auto-sampler	
Display options	Backlight	Yes	
	Auto Shutdown	1~60 min, off	
	IP Rating	IP54	
	Date and Time	Yes	
General	Electricity	AC Adapter, 100~240V AC input, DC24V output	
	Dimensions	220×195×68 mm	
	Weight	950g	

- BEP-M500T Multi-parameter Analyzer 1 Set**
E-301-QC 3 in 1 pH composite electrode 1 Piece
DJS-1VTC Conductivity electrode 1 Piece
DO-958-Q Polarographic DO electrode 1 Bottle
Standard pH Buffer Solution (50mL) (pH4.01, 7.00, 10.01) 1 Bottle each
Conductivity standard 1413µS/cm(50mL) 1 Bottle
REX-6 electrode stand 1 Set
GC65-240270-F Power adapter DC24V 1 Piece
Dust cover (Touch screen version) 1 Piece
European standard power cord(3VTJ2) 1 Piece
Universal USB cable 1 Piece
Touch pen (Capacitive screen) 1 Piece
Operation quick guide 1 Piece
Instruction manual 1 Piece



Benchtop Multiparameter Meter

BEP-M600L



Features

-  Color High contrast LCD touchscreen.
-  Data Storage 1000 sets for each parameter
-  IP54 waterproof.
-  Support for USB and RS-232 communication.



Description

- The BEP-M600L Benchtop Multi-parameter Analyzer adopts a multi-purpose measurement method, which is suitable for the fields of health disease control, environmental protection, biological agriculture and forestry, water quality analysis, scientific research and institutions of higher learning.

Features

General Features

- Color High contrast LCD touchscreen, 7.0 inches.
- Intelligent operation system provides features including user management, calibration management, method management, data management and etc.
- Multi-reading feature allows auto-read, timed-read and continuous-read.
- Automatic/Manual temperature compensation ensures accurate results.
- Data Storage 1000 sets for each parameter (GLP-compliant).
- Data analysis feature helps user review, compare and recalculate results.
- Support for USB and RS-232 communication.
- Reset feature automatically resumes all settings back to factory default options.
- IP54 waterproof.
- Support for three modules, namely pH/ISE module, conductivity module and dissolved oxygen module.
- Select from three modules to customize the four channels your way.
- Four channels can be configured in any combination.

pH

- 1-6 points calibration with Standard Recognition.
- Selectable pH buffer groups, including NIST, DIN, GB.
- Support for pH standard group customization.
- Automatic electrode diagnosis with pH slope and offset display.

Ion

- 1-6 points calibration
- Selectable measurement unit, including $\mu\text{g/L}$, mg/L , g/L , mmol/L , pX , etc.
- Multi-measurement modes are supported, including Direct Reading mode, Standard Addition mode, Sample Addition mode and GRAN mode.
- Over 10 methods are built-in, including F^- , Cl^- , Br^- , I^- , NO_3^- , BF_4^- , NH_4^+ , K^+ , Na^+ , Ca^{2+} , Cu^{2+} , Pb^{2+} , Ag^+ and etc., user-defined method is supported.

Conductivity

- 1-5 points calibration with Standard Recognition.
- Settable parameters, including cell constant, reference temperature (5/10/15/18/20/25°C), temperature compensation coefficient and TDS factor.
- Temperature compensation type (none, linear, pure water).

DO

- Support for air-saturated water or zero oxygen calibration.
- Auto barometric pressure compensation.
- Manual Salinity Factor Correction.
- Selectable pressure unit, including kPa, mbar, Torr, Atm.

Specification

Model			BEP-M600L
Parameters			pH/EC/ISE/DO/Temp. (mV/ORP/pX/Resistivity/TDS/Sal./DO Saturation)
pH	Range	-2.000 to 20.000 pH	
	Resolution	0.1, 0.01, 0.001pH	
	Accuracy	±0.002 pH	
	Calibration Points	Up to 6	
	Standard Customization	Yes	
	Calibration Reminder	Yes	
	Standard Recognition	NIST, GB, DIN buffers	
	Slope Limit	Yes	
mV/ORP	Range	-2000.00 to 2000.00 mV	
	Resolution	0.1, 0.01 mV	
	Accuracy	±0.03% FS	
	EH ORP Mode	Yes	
	ORP Calibration Points	1~2 custom point (relative mV)	
pX	Range	-2.000 to 20.000 pX	
	Resolution	0.1, 0.01, 0.001 pX	
	Accuracy	±0.002 pX	
	Calibration Points	Up to 6	
ISE	Range	1.000E-9 to 9.999E9	
	Units	mol/L, mmol/L, g/L, mg/L, µg/L	
	Resolution	Up to 4 significant digits	
	Accuracy	±0.3%	
	Calibration Points	Up to 6	
Conduc-tivity	Range	0.000 µS/cm to 3000 mS/cm	
	Resolution	0.001 µS/cm minimum, various with range selection	
	Accuracy	±0.5% FS	
	Reference Temperature	5, 10, 15, 18, 20, 25 °C	
	Calibration Points	Up to 5	
	Standard Recognition	146.5µS/cm, 1408µS/cm, 12.85mS/cm, 111.3mS/cm	
Resistivity	Range	5.00 Ω·cm~100.00 MΩ·cm	
	Resolution	0.01 Ω·cm minimum	
	Accuracy	±0.5% FS	
TDS	Range	0.000 mg/L~1000 g/L	
	Resolution	0.001ppm minimum, changed with range	
	Accuracy	±0.5%FS	
Salinity	Type	NaCl%	
	Range	(0.00~8.00)%	
	Resolution	0.01%	
	Accuracy	±0.1%	

Model			BEP-M600L
Dissolved	Range	0.00 to 99.99 mg/L	
Oxygen Concentration - Polarographic	Resolution	0.01 mg/L	
	Accuracy	±0.1 mg/L	
	Calibration Reminder	Yes	
	Calibration Points	Air-saturated water or zero point	
	Barometric Compensation	Yes	
	Manual Salinity Factor Correction	Yes	
% Saturation -Polarographic	Range	(0.0 to 600)%	
	Resolution	0.10%	
	Accuracy	±2.0%	
Temperature	Range	-5 to 130 °C	
	Unit	°C	
	Resolution	0.1	
	Accuracy	±0.1	
Measurement	Reading Mode	Auto Read (Fast, Medium, Slow), Timed, Continuous	
	Reading Prompts	Reading, Stable, Locked	
	Temp. Compensation	ATC, MTC	
Data management	Data Storage	1000 Groups	
	GLP Features	Yes	
	Log Management	Yes	
Inputs	pH Electrode	BNC(Q9)	
	DO with Temp. Probe	6-pin Mini-DIN	
	Conductivity with Temp. Probe	8-pin Mini-DIN	
Outputs	USB, RS 232	USB 2.0 flash memory device, printer, PC	
Display options	Backlight	Yes	
	Auto Shutdown	1~60 min, off	
	IP Rating	IP54	
	Date and Time	Yes	
General	Electricity	AC Adapter, 100~240V AC input, DC24V output	
	Dimensions	280×280×130 mm	
	Weight	2500g	

Description		Quantity	Description		Quantity
BEP-M600L Multi-parameter Analyzer		1 Set	GM65-240270 AC/DC Power adapter		1 pc
E-201-L pH composite electrode		1 pc	Touch pen		1 pc
T-818-L Temperature probe		1 pc	Dust cover		1 pc
PF-202-L Fluoride ion composite electrode		1 pc	Mini USB communication cable		1 pc
pH4.01/7.00/ 10.01Powder		3×5 sachets	Power cord		1 pc
DJS-1-L Conductivity electrode		1 pc	Ground wire(DC401)		1 pc
DO-958-L DO electrode		1 pc	Instruction manual		1 pc
REX-3 Electrode stand		1 Set	Operation guide		1 pc

Benchtop Multiparameter Meter

BEP-M900



Measurement Parameters

BEP-M900:
pH, mV, relative mV, ion concentration, conductivity, TDS, salinity, resistivity, DO, temperature.

Features

pH:
Multiparameter water quality meter is equipped with a 6.5 inches backlit LCD display
1 to 5 points calibration with automatic recognition for USA, NIST and DIN buffers
Automatic electrode diagnosis shows the pH slope and offset

ORP:
1 point offset calibration allows adjusting the displayed value to a known standard
Relative and absolute millivolt modes ensure the reliable ORP measurements

Ion Concentration:
2 to 5 points calibration, including the 8 concentration points can be selected
Direction concentration readout simplifies the measurement process
Automatic electrode diagnosis shows the calibration points and slopes
Selectable concentration units (ppm, mg/L, mol/L, mmol/L) and ionic valency



General Features



Auto-Read function senses and locks the measurement endpoint



Automatic temperature compensation ensures accurate readings over the entire range



Calibration due alarm prompts the user to calibrate the meter regularly



Setup menu allows setting the number of calibration points, resolution, stability criteria, temperature unit, date and time, etc.



Reset function automatically resumes all settings back to the factory defaults



Expanded memory stores or recalls up to 500 data sets
USB communication interface for data transfer and timed interval readings

Conductivity/TDS/Salinity/Resistivity:

1 to 5 points calibration with automatic recognition for conductivity standards

Selectable cell constant, reference temperature, TDS factor, linear and pure water compensations, seawater, and practical salinity measurement modes

Automatic electrode diagnosis shows the calibration points and factors



Dissolved Oxygen:

1 or 2 points calibration using the air-saturated water or zero oxygen solution

Salinity and barometric pressure compensations eliminate the measurement error



Ordering Information

BEP-M900:
pH/conductivity/dissolved oxygen electrodes, temperature probe, pH buffer sachets, conductivity standard solutions, DO electrolyte solution, membrane cap, electrode holder, USB cable, and power adapter.



Specifications

Model		BEP-M900
pH		
Range	-2.000~20.000pH	●
Resolution	0.001, 0.01, 0.1pH, selectable	●
Accuracy	±0.002pH	●
Calibration Points	1 to 5 points	●
pH Buffer Options	USA, NIST, DIN, 2 custom buffers	●
ORP		
Range	±1999.9mV	●
Resolution	0.1, 1mV, selectable	●
Accuracy	±0.2mV	●
Calibration Points	1 point	●
Ion		
Range	0.001~19999 (depending on the range of ISE)	●
Resolution	0.001, 0.01, 0.1, 1	●
Accuracy	±0.5% F.S. (monovalent), ±1% F.S. (divalent)	●
Measurement Units	ppm, mg/L, mol/L, mmol/L	●
Calibration Points	2 to 5 points (0.001, 0.01, 0.1, 1, 10, 100, 1000, 10000)	●
Conductivity		
Range	0.01~20.00, 200.0, 2000μS/cm, 20.00, 200.0mS/cm	●
Resolution	0.001,0.01,0.1,1	●
Accuracy	±0.5% F.S.	●
Calibration Points	1 to 5 points	●
Calibration Solutions	10μS/cm, 84μS/cm, 1413μS/cm, 12.88mS/cm, 111.8mS/cm	●
Temperature Coefficient	Linear (0.0~10.0%/°C), pure water	●
Reference Temperature	20/25°C	●
Cell Constant	K=0.1, 1, 10 or custom	●

Specifications

Model		BEP-M900
TDS		
Range	0~10.00, 100.0, 1000ppm, 10.00, 200.0ppt	●
Resolution	0.01,0.1,1	●
Accuracy	±1% F.S.	●
TDS Factor	0.1~1.0 (default 0.5)	●
Salinity		
Range	0.00~42.00psu, 0.00~80.00ppt	●
Resolution	0.01	●
Accuracy	±1% F.S.	●
Resistivity		
Range	0.00~20.00MΩ	●
Resolution	0.01	●
Accuracy	±1% F.S.	●
DO		
Range	0.00~20.00mg/L, 0.0~200.0% saturation	●
Resolution	0.01mg/L, 0.1%	●
Accuracy	±0.2mg/L, ±2.0%	●
Calibration Points	1or 2 points	●
Barometric Pressure Correction	60.0~112.5kPa/450~850mmHg, manual	●
Salinity Correction	0.0~50.0g/L, manual	●
General Spec.		
Temperature Compensation	0~100°C/32~212°F, manual or automatic	●
Memory	500 data sets, USB communication interface	●
Power Requirements	DC5V power adapter	●
Dimensions and Weight	210(L)×188(W)×60(H)mm,1.5kg	●

Portable Multiparameter Meter

BEP-M51X



- ◆ VA display screen for clearer visualization
- ◆ Lightweight, handheld design with simple button controls
- ◆ 5-minute auto-power-off function to extend battery life
- ◆ 5-in-1 Multifunctional Design—TDS, pH, EC, Salinity, Temperature

Applications

It is suitable for household tap water, outdoor water sources, aquariums, swimming pools, aquaculture, laboratory testing, laundry detergent residue, and water purifier filter replacement.

Features



Industrial-grade sensor solution ensures accurate readings with minimal error.



Hard cardboard box packaging with plastic tray for enhanced product protection.



5-minute auto-power-off function to extend battery life.



VA display screen for clearer visualization.



Acrylic screen panel for effective screen protection.

Specification

Model	BEP-M51X
TDS Range	0~9999ppm/10~200ppt
pH Range	0-14.00pH
EC Range	0~9990μS/cm/10~400mS/cm
Salinity Range	0~9999ppm/10~200ppt
Temperature Range	0-60 °C /32-140 °F
Resolution	1ppm/0.01pH/1μS/cm/1ppm
Accuracy	0~999ppm (±2%);1000~9990ppm (±4%)
Battery Type	LR44
Battery Quantity	4pcs.
Operating Voltage	6V
Operating Current	12mA
Operating Environment	0°C~60°C
IP Rating	IP65
Shell Material	ABS
External Dimensions	185*40*40mm
Package Dimensions	216.5*65*43.5mm
N.W./G.W.	110g/180g
Quantity per Carton	58pcs.
Carton Gross Weight	10.3kg

Portable Multiparameter Meter

BEP-M71X



- ◆ VA display screen for clearer visualization
- ◆ Lightweight, handheld design with simple button controls
- ◆ 5-minute auto-power-off function to extend battery life
- ◆ 7-in-1 Multifunctional Design—TDS, pH, EC, Salinity, S.G, ORP, Temperature

Applications

It is suitable for household tap water, outdoor water sources, aquariums, swimming pools, aquaculture, laboratory testing, laundry detergent residue, and water purifier filter replacement.

Features



Industrial-grade sensor solution ensures accurate readings with minimal error.



Hard cardboard box packaging with plastic tray for enhanced product protection.



5-minute auto-power-off function to extend battery life.



VA display screen for clearer visualization.



Acrylic screen panel for effective screen protection.

Specification

Model	BEP-M71X
TDS Range	0~9999ppm/10~200ppt
pH Range	0~14.00pH
EC Range	0~9990μS/cm/10~400mS/cm
Salinity Range	0~9999ppm/10~200ppt
S.G Range	0~1.222
ORP Range	±990mV
Temperature Range	0~60℃/32~140℉
Resolution	1 ppm/ 0.01pH/ 1μS/cm/ 1ppm/ 1mV
Accuracy	0~999ppm (±2%); 1000~9990ppm (±4%)
Battery Type	LR44
Battery Quantity	4pcs.
Operating Voltage	6V
Operating Current	12mA
Operating Environment	0℃~60℃
IP Rating	IP65
Shell Material	ABS
External Dimensions	185*40*40mm
Package Dimensions	216.5*65*43.5mm
N.W./G.W.	110g/180g
Quantity per Carton	58pcs.
Carton Gross Weight	10.3kg

Portable Multiparameter Meter

BEP-M310A



pH

- ◆ 1-5 points calibration with Standard Recognition.
- ◆ Selectable pH buffer groups, including NIST, DIN, GB.
- ◆ Automatic electrode diagnosis with pH slope and offset display.

Conductivity

- ◆ 1-3 points calibration with Standard Recognition.
- ◆ Settable parameters, including cell constant, temperature compensation coefficient, and TDS factor.
- ◆ Temperature compensation type (none, linear, pure water).

Ion

- ◆ 1-5 points calibration.
- ◆ Selectable measurement unit, including $\mu\text{g/L}$, mg/L , g/L , mmol/L , pX , ppm , etc.
- ◆ Measurement modes are supported, including Direct Reading mode, Standard Addition mode, Sample Addition mode and GRAN mode.
- ◆ Over 10 methods are built-in, including F^- , Cl^- , Br^- , I^- , NO_3^- , BF_4^- , NH_4^+ , K^+ , Na^+ , Ca^{2+} , Cu^{2+} , Pb^{2+} , Ag^+ and etc., user-defined method is supported.

Features



LCD display screen, 3.5 inches.



Support for USB communication. Multiple external devices available
PC, printer



Tilting stand
Support better reading and operation



Wristbands Safety hold

Description

The BEP-M310A Portable Multi-parameter Analyzer adopts a multi-purpose measurement method, which is suitable for the fields of health disease control, environmental protection, biological agriculture and forestry, water quality analysis, scientific research and institutions of higher learning.

General Features



IP65 waterproof.



Data capacity of up to 500 sets for each parameter (GLP-compliant).



Multi-reading feature allows auto-read, timed-read and continuous-read.



Auto-power off feature effectively extends the battery service life.
Reset feature automatically resumes all settings back to factory default options.



Automatic/Manual temperature compensation ensures accurate results.
Auto-hold feature senses and locks the measurement endpoint.

Specifications

Model		BEP-M310A
Parameters		pH/EC/ISE/Temp. (mV/ORP/pX/Resistivity/TDS/Sal.)
pH	Range	-2.00 to 20.00 pH
	Resolution	0.1, 0.01 pH
	Accuracy	±0.01 pH
	Calibration Points	Up to 5
	Standard Customization	Yes
	Standard Recognition	NIST, GB and DIN buffers
	Slope Limit	Yes
mV	Range	-2000.0 to 2000.0 mV
	Resolution	0.1
	Accuracy	±0.3 mV or ±0.1%
pX	Range	-2.00 to 20.00
	Resolution	0.1, 0.01 pX
	Accuracy	±0.01 pX
	Calibration Points	Up to 5
ISE	Range	1E-9 to 9.999E9
	Units	mol/L, mmol/L, g/L, mg/L, µg/L
	Resolution	Up to 4 significant digits
	Accuracy	±0.5%
	Calibration Points	Up to 5
Conductivity	Range	0.000 µS/cm to 500 mS/cm
	Resolution	0.001 µS/cm minimum,various with range selection
	Accuracy	±1.0% FS
	Reference Temperature	20, 25 °C
	Calibration Points	Up to 3
	Standard Recognition	84µS/cm, 1413µS/cm, 12.88mS/cm
Resistivity	Range	5.00 Ω·cm~20.00 MΩ·cm
	Resolution	0.01 Ω·cm minimum
	Accuracy	±1.0% FS
TDS	Range	0.00 ppm~300 ppt
	Resolution	0.01ppm minimum, various with range selection
	Accuracy	±1.0%FS
Salinity	Range	0.0~80.0 ppt
	Resolution	0.1ppt
	Accuracy	±2ppt

Model		BEP-M310A
Temperature	Range	-5 to 110 °C, 23 to 230 °F
	Unit	°C, °F
	Resolution	0.1
	Relative Accuracy	±0.2
Measurement	Reading Mode	Auto Read (Fast, Medium, Slow), Timed, Continuous
	Reading Prompts	Reading, Stable, Locked
	Temp. Compensation	ATC, MTC
Data manage-ment	Data Storage	500 results each
	GLP Features	Yes
Inputs	pH Electrode	BNC(Q9)
Inputs	Conductivity with Temp. Probe	5-pin aviation connector
Outputs	USB	PC, printer
Display options	Backlight	Yes
	Auto Shutdown	300, 600, 1200, 1800, 3600 sec., off
	IP Rating	IP65
	Date and Time	Yes
General	Electricity	Rechargeable Lithium battery, AC Adapter, 100-240V AC input, DC5V output
	Dimensions	W80×D255×H35 mm
	Package Dimension (W*D*H)(mm)	W360×D300×H180
	N.W./G.W.(kg)	0.4kg/2.8kg

BEP-M310A Multi-parameter Analyzer 1 Set

- MicroUSB to DB25 printing cable 1 Piece
- Silicone rubber case (Installed on the host) 1 Piece
- Wristband assembly(Small) 1 Piece
- Carrying case(Large) 1 Piece
- Operation quick guide 1 Piece
- Instruction manual 1 Piece



- E-301-QC 3 in 1 pH composite electrode 1 Piece
- DJS-IVTC Conductivity electrode 1P iece
- Standard pH Buffer Solution (50mL) (pH4.01, 7.00, 10.01) 1 Bottle each
- Conductivity standards 1413µS/cm(50mL) 1 Bottle
- Power adapter DC 5V 2A(European plug model, GAT-0502000) 1 Piece
- Electrode holder 2 Pieces
- MicroUSB cable 1 Piece



Portable Multiparameter Meter

BEP-M310F



The BEP-M310F Portable Multi-parameter Analyzer adopts a multi-purpose measurement method, which is suitable for the fields of health disease control, environmental protection, biological agriculture and forestry, water quality analysis, scientific research and institutions of higher learning.

Features

- LCD display screen, 3.5 inches.
- Multi-reading feature allows auto-read, timed-read and continuous-read.
- Automatic/Manual temperature compensation ensures accurate results.
- Auto-hold feature senses and locks the measurement endpoint.
- Data capacity of up to 500 sets for each parameter (GLP-compliant).
- Support for USB communication.
- Auto-power off feature effectively extends the battery service life.
- Reset feature automatically resumes all settings back to factory default options.
- IP65 waterproof.

Features

- ◆ **pH**
1-5 points calibration with Standard Recognition.
Selectable pH buffer groups, including NIST, DIN, GB.
Automatic electrode diagnosis with pH slope and offset display.
- ◆ **Ion**
1-5 points calibration.
Selectable measurement unit, including $\mu\text{g/L}$, mg/L , g/L , mmol/L , pX , ppm , etc.
Measurement modes are supported, including Direct Reading mode, Standard Addition mode, Sample Addition mode and GRAN mode.
Over 10 methods are built-in, including F^- , Cl^- , Br^- , I^- , NO_3^- , BF_4^- , NH_4^+ , K^+ , Na^+ , Ca^{2+} , Cu^{2+} , Pb^{2+} , Ag^+ and etc., user-defined method is supported.
- ◆ **Conductivity**
1-3 points calibration with Standard Recognition.
Settable parameters, including cell constant, temperature compensation coefficient, and TDS factor.
Temperature compensation type (none, linear, pure water).
- ◆ **DO**
Support for air-saturated water or zero oxygen calibration.
Auto barometric pressure compensation
Manual Salinity Factor Correction
Selectable pressure unit, including kPa , mbar , Torr , Atm .



LCD display screen.
Data capacity of up to 500 sets for each parameter.
Support for USB communication.
IP65 waterproof.

Specifications

Model		BEP-M310F
Parameters		pH/EC/ISE/DO/Temp. (mV/ORP/pX/Resistivity/TDS/Sal./DO Saturation)
pH	Range	-2.00 to 20.00pH
	Resolution	0.1, 0.01 pH
	Accuracy	±0.01 pH
	Calibration Points	Up to 5
	Standard Customization	Yes
	Standard Recognition	NIST, GB and DIN buffers
	Slope Limit	Yes
mV	Range	-2000.0 to 2000.0 mV
	Resolution	0.1
	Accuracy	±0.3 mV or ±0.1%
pX	Range	-2.00 to 20.00
	Resolution	0.1, 0.01 pX
	Accuracy	±0.01 pX
	Calibration Points	Up to 5
ISE	Range	1E-9 to 9.999E9
	Units	mol/L, mmol/L, g/L, mg/L, µg/L
	Resolution	Up to 4 significant digits
	Accuracy	±0.5%
	Calibration Points	Up to 5
Conductivity	Range	0.000 µS/cm to 500 mS/cm
	Resolution	0.001 µS/cm minimum,various with range selection
	Accuracy	±1.0% FS
	Reference Temperature	20, 25 °C
	Calibration Points	Up to 3
	Standard Recognition	84µS/cm, 1413µS/cm, 12.88mS/cm
Resistivity	Range	5.00 Ω·cm~20.00 MΩ·cm
	Resolution	0.01 Ω·cm minimum
	Accuracy	±1.0% FS
TDS	Range	0.00 ppm~300 ppt
	Resolution	0.01ppm minimum, various with range selection
	Accuracy	±1.0%FS

BEP-M310F Multi-parameter Analyzer 1 Set

- Electrode holder 2 Pieces
- MicroUSB cable 1 Piece
- MicroUSB to DB25 printing cable 1 Piece
- Silicone rubber case (Installed on the host) 1 Piece
- Wristband assembly(Small) 1 Piece
- Carrying case(Large) 1 Piece
- Operation quick guide 1 Piece
- Instruction manual 1 Piece

- E-301-QC 3 in 1 pH composite electrode 1 Piece
- DJS-1VTC Conductivity electrode 1 Piece
- DO-958-Q Polarographic DO electrode 1 Piece
- Standard pH Buffer Solution (50mL) (pH4.01, 7.00, 10.01) 1 Bottle each
- Conductivity standards 1413µS/cm(50mL) 1 Bottle
- Power adapter DC 5V 2A(European plug model, GAT-0502000) 1 Piece

Specifications

Model		BEP-M310F
Salinity	Range	0.0~80.0 ppt
	Resolution	0.1ppt
	Accuracy	±2ppt
Dissolved Oxygen Con- centration	Range	0.00 to 50.00ppm
	Resolution	0.01ppm
	Accuracy	±0.1ppm
	Calibration Points	Air-saturated water or zero point
	Barometric Compensation	Yes
	Manual Salinity Factor Correction	Yes
Saturation	Range	(0.0 to 300.0)%
	Resolution	0.1%
	Accuracy	±10.0%
Temperature	Range	-5 to 110 °C, 23 to 230 °F
	Unit	°C, °F
	Resolution	0.1
	Relative Accuracy	±0.2
Measurement	Reading Mode	Auto Read (Fast, Medium, Slow), Timed, Continuous
	Reading Prompts	Reading, Stable, Locked
	Temp. Compensation	ATC, MTC
Data manage- ment	Data Storage	500 results each
	GLP Features	Yes
Inputs	pH Electrode	BNC(Q9)
	DO with Temp. Probe	4-pin aviation connector
	Conductivity with Temp. Probe	5-pin aviation connector
Outputs	USB	PC, printer
Display options	Backlight	Yes
	Auto Shutdown	300, 600, 1200, 1800, 3600 sec., off
	IP Rating	IP65
	Date and Time	Yes
General	Electricity	Rechargeable Lithium battery, AC Adapter, 100-240V AC input, DC5V output
	Dimensions	80×255×35 mm
	Weight	400g



Portable Multiparameter Meter

BEP-M310T



Features

- pH**
1-5 points calibration with Standard Recognition.
Selectable pH buffer groups, including NIST, DIN, GB.
Automatic electrode diagnosis with pH slope and offset display.
- Ion**
1-5 points calibration.
Selectable measurement unit, including $\mu\text{g/L}$, mg/L , g/L , mmol/L , pX , ppm , etc.
Measurement modes are supported, including Direct Reading mode, Standard Addition mode, Sample Addition mode and GRAN mode.
Over 10 methods are built-in, including F^- , Cl^- , Br^- , I^- , NO_3^- , BF_4^- , NH_4^+ , K^+ , Na^+ , Ca^{2+} , Cu^{2+} , Pb^{2+} , Ag^+ and etc., user-defined method is supported.
- Conductivity**
1-3 points calibration with Standard Recognition.
Settable parameters, including cell constant, temperature compensation coefficient, and TDS factor.
Temperature compensation type (none, linear, pure water).
- DO**
Support for air-saturated water or zero oxygen calibration.
Auto barometric pressure compensation
Manual Salinity Factor Correction
Selectable pressure unit, including kPa , mbar , Torr , Atm .



The BEP-M310T Portable Multi-parameter Analyzer adopts a multi-purpose measurement method, which is suitable for the fields of health disease control, environmental protection, biological agriculture and forestry, water quality analysis, scientific research and institutions of higher learning.

Features

- 1.1LCD display screen,4.3 inches.
- Multi-reading feature allows auto-read, timed-read and continuous-read.
- Automatic/Manual temperature compensation ensures accurate results.
- Auto-hold feature senses and locks the measurement endpoint.
- Data capacity of up to 500 sets for each parameter (GLP-compliant).
- Support for USB communication.
- Auto-power off feature effectively extends the battery service life.
- Reset feature automatically resumes all settings back to factory default options.
- IP65 waterproof.



LCD display screen.
Data capacity of up to 500 sets for each parameter.
Support for USB communication.
IP65 waterproof.

Specifications

Model		BEP-M310T
Parameters		pH/EC/ISE/DO/Temp. (mV/pX/Resistivity/TDS/Sal./DO Saturation)
pH	Range	-2.00 to 20.00pH
	Resolution	0.1, 0.01 pH
	Accuracy	±0.01 pH
	Calibration Points	Up to 5
	Standard Customization	Yes
	Calibration Reminder	Yes
	Standard Recognition	NIST, GB and DIN buffers
	Slope Limit	Yes
mV	Range	-2000.0 to 2000.0 mV
	Resolution	0.1 mV
	Accuracy	±0.3 mV or ±0.1%
pX	Range	-2.00 to 20.00 pX
	Resolution	0.1, 0.01 pX
	Accuracy	±0.01 pX(±0.01 pX)
	Calibration Points	Up to 5
ISE	Range	1E-9 to 9.999E9
	Units	mol/L, mmol/L, g/L, mg/L, µg/L
	Resolution	Up to 4 significant digits
	Accuracy	±0.5%
	Calibration Points	Up to 5
Conductivity	Range	0.000 µS/cm to 1000 mS/cm
	Resolution	0.001 µS/cm minimum,various with range selection
	Accuracy	±1.0% FS
	Reference Temperature	20, 25 °C
	Calibration Points	Up to 3
	Standard Recognition	84µS/cm, 1413µS/cm, 12.88mS/cm
Resistivity	Range	5.00 Ω ·cm~20.00 MΩ ·cm
	Resolution	0.01Ω ·cm minimum,
	Accuracy	±1.0% FS
TDS	Range	0.00ppm~300ppt
	Resolution	0.01mg/L minimum, 3 significant digits
	Accuracy	±1.0%

BEP-M310T Multi-parameter Analyzer 1 Set MicroUSB cable 1 Piece MicroUSB to DB25 printing cable 1 Piece Touch pen (Capacitive screen) 1 Piece Silicone rubber case (Installed on the host) 1 Piece Wristband assembly(Large) 1 Piece Carrying case(Large) 1 Piece Operation quick guide 1 Piece Instruction manual 1 Piece	E-301-QC 3 in 1 pH composite electrode 1 Piece DJS-IVTC Conductivity electrode 1 Piece DO-958-Q Polarographic DO electrode 1 Piece Standard pH Buffer Solution (50mL) (pH4.01, 7.00, 10.01) 1 Bottle each Conductivity standards 1413µS/cm(50mL) 1 Bottle Power adapter DC 5V 2A(European plug model, GAT-0502000) 1 Piece Electrode holder 2 Pieces
---	--

Model		BEP-M310T
Salinity	Range	0.0~80.0ppt
	Resolution	0.1ppt
	Accuracy	±2ppt
Dissolved Oxygen Con- centration	Sensor Type	Polarographic
	Range	0.00 to 50.00 ppm
	Resolution	0.01ppm
	Accuracy	±0.1ppm
	Calibration Points	Air-saturated water or zero point
	Barometric Compensation	Yes
	Manual Salinity Factor Correction	Yes
Saturation	Range	(0.0 to 300.0)%
	Resolution	0.10%
	Accuracy	±2.0%
Temperature	Range	-5 to 110 °C, 23.0 to 230 °F
	Unit	°C, °F
	Resolution	0.1
	Relative Accuracy	±0.2
Measurement	Reading Mode	Auto-Read (Fast, Medium, Slow), Timed, Continuous
	Reading Prompts	Reading, Stable, Locked
	Temp. Compensation	ATC, MTC
Data manage- ment	Data Storage	1000 Groups
	GLP Features	Yes
Inputs	pH Electrode	BNC(Q9)
	DO with Temp. Probe	4-pin aviation connector
	Conductivity with Temp. Probe	5-pin aviation connector
Outputs	USB	PC, printer
	Backlight	Yes
Display options	Auto Shutdown	300, 600, 1200, 1800, 3600sec., off
	IP Rating	IP65
	Date and Time	Yes
General	Electricity	Rechargeable Lithium battery, AC Adapter, 100-240V AC input, DC5V output
	Dimensions	90 x 255 x 40 mm
	Weight	500g



Portable Multiparameter Meter

BEP-M510A



pH

- ◆ 1-8 points calibration with Standard Recognition.
- ◆ Selectable pH buffer groups, including USA, NIST, DIN, GB, MERK, JIS.
- ◆ User-defined pH buffer or buffer group is supported.
- ◆ Automatic electrode diagnosis with pH slope and offset display.

Conductivity

- ◆ 1-5 points calibration automatically recognizes the standard solutions.
- ◆ Settable parameters, including cell constant, reference temperature (5/10/15/18/20/25°C), temperature compensation coefficient and TDS factor.
- ◆ Temperature compensation type (none, linear, pure water).

Ion

- ◆ 1-6 points calibration.
- ◆ Selectable measurement unit, including $\mu\text{g/L}$, mg/L , g/L , mmol/L , pX , ppm , etc.
- ◆ Multi-measurement modes are supported, including Direct Reading mode, Standard Addition mode, Sample Addition mode and GRAN mode.
- ◆ Over 10 methods are built-in, including F^- , Cl^- , Br^- , I^- , NO_3^- , BF_4^- , NH_4^+ , K^+ , Na^+ , Ca^{2+} , Cu^{2+} , Pb^{2+} , Ag^+ and etc., user-defined method is supported.

Features



Color High contrast LCD touchscreen, 4.3 inches.



Support for USB communication. Multiple external devices available
PC, printer



Tilting stand
Support better reading and operation



Wristbands Safety hold

Description

The BEP-M510A Portable Multi-parameter Analyzer adopts a multi-purpose measurement method, which is suitable for the fields of health disease control, environmental protection, biological agriculture and forestry, water quality analysis, scientific research and institutions of higher learning.

General Features



IP65 waterproof.



Intelligent operation system provides features including user management, calibration management, electrode management, method management, data management, Log Management and etc.



Data Storage 1000 sets for each parameter (GLP-compliant).
Data analysis feature helps user review, compare and recalculate results.



Multi-reading feature allows auto-read, timed-read and continuous-read.



Automatic/Manual temperature compensation ensures accurate results.
Auto-hold feature senses and locks the measurement endpoint.



Auto-power off feature effectively extends the battery service life.
Reset feature automatically resumes all settings back to factory default options.

Specifications

Model		BEP-M510A
Parameters		pH/EC/ISE/Temp. (mV/ORP/pX/Resistivity/TDS/Sal.)
pH	Range	-2.000 to 20.000 pH
	Resolution	0.1, 0.01, 0.001pH
	Accuracy	±0.002 pH
	Calibration Points	Up to 8
	Standard Customization	Yes
	Calibration Reminder	Yes
	Standard Recognition	USA, NIST, MERK, JIS, GB and DIN buffers
	Slope Limit	Yes
mV/ORP	Range	-2000.00 to 2000.00 mV
	Resolution	0.1, 0.01 mV
	Accuracy	±0.1 mV or ±0.03%
	EH ORP Mode	Yes
	ORP Calibration Points	1 custom point (relative mV)
pX	Range	-2.000 to 20.000 pX
	Resolution	0.1, 0.01, 0.001 pX
	Accuracy	±0.002 pX
	Calibration Points	Up to 6
ISE	Range	1E-9 to 9.999E9
	Units	mol/L, mmol/L, g/L, mg/L, µg/L
	Resolution	Up to 4 significant digits
	Accuracy	±0.3%
	Calibration Points	Up to 8
Conductivity	Range	0.000 µS/cm to 3000 mS/cm
	Resolution	0.001 µS/cm minimum,various with range selection
	Accuracy	±0.5% FS
	Reference Temperature	5, 10, 15, 18, 20, 25 °C
	Calibration Points	Up to 5
	Calibration Reminder	Yes
	Standard Recognition	10µS/cm, 84µS/cm, 500µS/cm, 1413µS/cm, 12.88mS/cm
Resistivity	Range	5.00 Ω·cm~100.00 MΩ·cm
	Resolution	0.01 Ω·cm minimum
	Accuracy	±0.5% FS
TDS	Range	0.000ppm~1000ppt
	Resolution	0.001ppm minimum, changed with range
	Accuracy	±0.5%FS

Model		BEP-M510A
Salinity	Range	0.0~80.0ppt
	Resolution	0.1ppt
	Accuracy	±1ppt
Temperature	Range	-10 to 135 °C, 14 to 275 °F
	Unit	°C, °F
	Resolution	0.1
	Relative Accuracy	±0.1
Measurement	Reading Mode	Auto Read (Fast, Medium, Slow), Timed, Continuous
	Reading Prompts	Reading, Stable, Locked
	Temp. Compensation	ATC, MTC
Data manage- ment	Data Storage	1000 Groups
	GLP Features	Yes
Inputs	pH Electrode	BNC(Q9)
	DO with Temp. Probe	4-pin aviation connector
	Conductivity with Temp. Probe	5-pin aviation connector
Outputs	USB	USB 2.0 flash memory device, PC, scanner
	Backlight	Yes
Display options	Auto Shutdown	300, 600, 1200, 1800, 3600 sec., off
	IP Rating	IP65
	Date and Time	Yes
General	Electricity	Rechargeable Lithium battery, AC Adapter, 100~240V AC input, DC5V output
	Dimensions	90×255×40 mm
	Package Dimension (W*D*H)	360×300×180mm
	G.W./N.W.	2.6kg/0.5kg

BEP-M510A Multi-parameter Analyzer 1 Set

OTG adapter cable (MicroUSB connector) 1 Piece
Touch pen (Capacitive screen) 1 Piece
Silicone rubber case (Installed on the host) 1 Piece
Wristband assembly(Large) 1 Piece
Carrying case(Large) 1 Piece
Operation quick guide 1 Piece
Instruction manual 1 Piece



E-301-QC 3 in 1 pH composite electrode 1 Piece
DJS-IVTC Conductivity electrode 1 Piece
Standard pH Buffer Solution (50mL) (pH4.01, 7.00, 10.01) 1 Bottle each
Conductivity standard 1413µS/cm(50mL) 1 Bottle
Power adapter DC 5V 2A(European plug model, GAT-0502000) 1 Piece
Electrode holder 2 Pieces
MicroUSB cable 1 Piece



Portable Multiparameter Meter

BEP-M510T




The BEP-M510T Portable Multi-parameter Analyzer adopts a multi-purpose measurement method, which is suitable for the fields of health disease control, environmental protection, biological agriculture and forestry, water quality analysis, scientific research and institutions of higher learning.

Features

- Color High contrast LCD touchscreen, 4.3 inches.
- Intelligent operation system provides features including user management, calibration management, electrode management, method management, data management, Log Management and etc.
- Multi-reading feature allows auto-read, timed-read and continuous-read.
- Automatic/Manual temperature compensation ensures accurate results.
- Auto-hold feature senses and locks the measurement endpoint.
- Data Storage 1000 sets for each parameter (GLP-compliant).
- Data analysis feature helps user review, compare and recalculate results.
- Support for USB communication.
- Auto-power off feature effectively extends the battery service life.
- Reset feature automatically resumes all settings back to factory default options.
- IP65 waterproof.

Features

- **pH**
1-8 points calibration with Standard Recognition.
Selectable pH buffer groups, including USA, NIST, DIN, GB, MERK, JIS.
User-defined pH buffer or buffer group is supported.
Automatic electrode diagnosis with pH slope and offset display.
- **Ion**
1-6 points calibration.
Selectable measurement unit, including $\mu\text{g/L}$, mg/L , g/L , mmol/L , pX , ppm , etc.
Multi-measurement modes are supported, including Direct Reading mode, Standard Addition mode, Sample Addition mode and GRAN mode.
Over 10 methods are built-in, including F^- , Cl^- , Br^- , I^- , NO_3^- , BF_4^- , NH_4^+ , K^+ , Na^+ , Ca^{2+} , Cu^{2+} , Pb^{2+} , Ag^+ and etc., user-defined method is supported.
- **Conductivity**
1-5 points calibration automatically recognizes the standard solutions.
Settable parameters, including cell constant, reference temperature ($5/10/15/18/20/25^\circ\text{C}$), temperature compensation coefficient and TDS factor.
Temperature compensation type (none, linear, pure water).
- **DO**
Support for air-saturated water or zero oxygen calibration.
Auto barometric pressure compensation
Manual Salinity Factor Correction
Selectable pressure unit, including kPa , mbar , Torr , Atm .

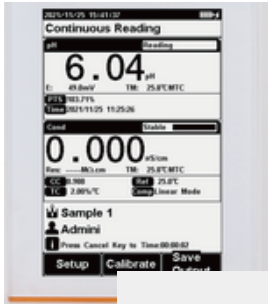

Color High contrast LCD touchscreen.
Data Storage 1000 sets for each parameter.
Support for USB communication.
IP65 waterproof.

Specifications

Model		BEP-M510T
Parameters		pH/EC/ISE/DO/Temp. (mV/ORP/pX/Resistivity/TDS/Sal./DO Saturation)
pH	Range	-2.000 to 20.000 pH
	Resolution	0.1, 0.01, 0.001pH
	Accuracy	±0.002 pH
	Calibration Points	Up to 8
	Standard Customization	Yes
	Calibration Reminder	Yes
	Standard Recognition	USA, NIST, MERK, JIS, GB and DIN buffers
	Slope Limit	Yes
mV/ORP	Range	-2000.00 to 2000.00 mV
	Resolution	0.1, 0.01 mV
	Accuracy	±0.1 mV or ±0.03%
	EH ORP Mode	Yes
pX	ORP Calibration Points	1 custom point (relative mV)
	Range	-2.000 to 20.000 pX
	Resolution	0.1, 0.01, 0.001 pX
	Accuracy	±0.002 pX
Calibration Points	Up to 6	
ISE	Range	1E-9 to 9.999E9
	Units	mol/L, mmol/L, g/L, mg/L, µg/L
	Resolution	Up to 4 significant digits
	Accuracy	±0.3%
Calibration Points	Up to 8	
Conductivity	Range	0.000 µS/cm to 3000 mS/cm
	Resolution	0.001 µS/cm minimum,various with range selection
	Accuracy	±0.5% FS
	Reference Temperature	5, 10, 15, 18, 20, 25 °C
	Calibration Points	Up to 5
	Calibration Reminder	Yes
Standard Recognition	10µS/cm, 84µS/cm, 500µS/cm, 1413µS/cm, 12.88mS/cm	
Resistivity	Range	5.00 Ω·cm~100.00 MΩ·cm
	Resolution	0.01 Ω·cm minimum
	Accuracy	±0.5% FS
TDS	Range	0.000ppm~1000ppt
	Resolution	0.001ppm minimum, changed with range
	Accuracy	±0.5%FS

- BEP-M510T Multi-parameter Analyzer 1 Set**
MicroUSB cable 1 Piece
OTG adapter cable (MicroUSB connector) 1 Piece
Touch pen (Capacitive screen) 1 Piece
Silicone rubber case (Installed on the host) 1 Piece
Wristband assembly(Large) 1 Piece
Carrying case(Large) 1 Piece
Operation quick guide 1 Piece
- Instruction manual 1 Piece
E-301-QC 3 in 1 pH composite electrode 1 Piece
DJS-IVTC Conductivity electrode 1 Piece
DO-958-Q Polarographic DO electrode 1 Piece
Standard pH Buffer Solution (50mL) (pH4.01, 7.00, 10.01) 1 Bottle each
Conductivity standard 1413µS/cm(50mL) 1 Bottle
Power adapter DC 5V 2A(European plug model, GAT-0502000) 1 Piece
Electrode holder 2 Pieces

Model		BEP-M510T
Salinity	Range	0.0~80.0ppt
	Resolution	0.1ppt
	Accuracy	±1ppt
Dissolved Oxygen Con- centration	Sensor Type	Polarographic
	Range	0.00 to 99.99 ppm
	Resolution	0.01ppm
	Accuracy	±0.1ppm (0.00 to 50.00 ppm);±0.5ppm (50.00 to 99.99 ppm)
	Calibration Reminder	Yes
	Calibration Points	Air-saturated water or zero point
	Barometric Compensation	Yes
	Manual Salinity Factor Correction	Yes
Saturation	Range	(0.0 to 600.0)%
	Resolution	0.10%
	Accuracy	±10.0%
Temperature	Range	-10 to 135 °C, 14 to 275 °F
	Unit	°C, °F
	Resolution	0.1
	Relative Accuracy	±0.1
Measurement	Reading Mode	Auto Read (Fast, Medium, Slow), Timed, Continuous
	Reading Prompts	Reading, Stable, Locked
	Temp. Compensation	ATC, MTC
Data manage- ment	Data Storage	1000 Groups
	GLP Features	Yes
Inputs	pH Electrode	BNC(Q9)
	DO with Temp. Probe	4-pin aviation connector
	Conductivity with Temp. Probe	5-pin aviation connector
Outputs	USB	USB 2.0 flash memory device, PC, scanner
Display options	Backlight	Yes
	Auto Shutdown	300, 600, 1200, 1800, 3600 sec., off
	IP Rating	IP65
	Date and Time	Yes
General	Electricity	Rechargeable Lithium battery, AC Adapter, 100~240V AC input, DC5V output
	Dimensions	90×255×40 mm
	Weight	500g



Portable Multiparameter Meter

BEP-M900P



Measurement Parameters

BEP-M900P: pH, mV, relative mV, ion concentration, conductivity, TDS,salinity, resistivity, DO, temperature.

General Features

Calibration due alarm prompts user to calibrate the meter regularly

Expanded memory stores or recalls up to 500 data sets

USB communication interface is easy to transfer data to PC

Multi-mode power scheme (batteries, power adapter and USB port) ensures that using

Setup menu allows setting the number of calibration points, resolution, stability criteria, temperature unit, date and time, etc.

Reset function automatically resumes all settings back to the factory defaults

Automatic temperature compensation ensures accurate readings over the entire range
Auto-Read function senses and locks the measurement endpoint
Auto-Power Off effectively conserves battery life

Features

pH:
Multi parameter water quality meter is equipped with a 3.5 inches backlit LCD display
1 to 5 points calibration with automatic recognition for USA, NIST and DIN buffers
Automatic electrode diagnosis shows the pH slope and offset

ORP:
1 point offset calibration allows adjusting the displayed value to a known standard
Relative and absolute millivolt modes ensure the reliable ORP measurements

Dissolved Oxygen:
1 or 2 points calibration using the air-saturated water or zero oxygen solution
Salinity and barometric pressure compensations eliminate the measurement error

Ion Concentration:
2 to 5 points calibration, including the 8 concentration points can be selected
Direct ion concentration readout simplifies the measurement process
Automatic electrode diagnosis shows the calibration points and slopes
Selectable concentration units (ppm, mg/L, mol/L, mmol/L) and ionic valency

Conductivity/TDS/Salinity/Resistivity:
1 to 5 points calibration with automatic recognition for conductivity standards
Selectable cell constant, reference temperature, TDS factor, linear and pure water compensations, seawater and practical salinity measurement modes
Automatic electrode diagnosis shows the calibration points and factors



Specifications

Model		BEP-M900P
pH		
Range	-2.000~20.000pH	●
Resolution	0001,0.01,0.1pH, selectable	●
Accuracy	0.002pH	●
Calibration Points	1 to 5 points	●
pHBuffer Options	USA, NIST, DIN, 2 custom buffers	●
ORP		
Range	±1999.9mV	●
Resolution	0.1, 1mV, selectable	●
Accuracy	±0.2mV	●
Calibration Points	1 point	●
Ion		
Range	0.001~19999 (depending on the range of ISE)	●
Resolution	0.001, 0.01, 0.1, 1	●
Accuracy	±0.5% F.S. (monovalent), ±1% F.S. (divalent)	●
Measurement Units	ppm, mg/L, mol/L, mmol/L	●
Calibration Paints	2 to 5 points (0.001, 0.01, 0.1, 1, 10, 100, 1000, 10000)	●
Conductivity		
Range	0.01~20.00, 200.0, 2000μS/cm, 20.00, 200.0mS/cm	●
Resolution	0001,0.01,0.1,1	●
Accuracy	±0.5% F.S.	●
Calibration Paints	1 to 5 points	●
Calibration Solutions	10μS/cm, 84μS/cm, 1413μS/cm, 12.88mS/cm, 111.8mS/cm	●
Temperature Coefficient	Linear(0.0~10.0%/°C).pure water	●
Reference Temperature	20/25°C	●
Cell Constant	K=0.1, 1, 10 or custom	●

Specifications

Model		BEP-M900P
TDS		
Range	0~10.00, 100.0, 1000ppm, 10.00, 200.0ppt	●
Resolution	001,0.1,1	●
Accuracy	±1% F.S.	●
TDS Factor	0.1~1.0 (default 0.5)	●
Salinity		
Range	0.00~42.00psu, 0.00~80.00ppt	●
Resolution	0.01	●
Accuracy	±1% F.S.	●
Resistivity		
Range	0.00~20.00MΩ	●
Resolution	0.01	●
Accuracy	±1% F.S.	●
DO		
Range	0.00~20.00mg/L, 0.0~200.0% saturation	●
Resolution	001mg/L,0.1%	●
Accuracy	±0.2mg/L, ±2.0%	●
Calibration Points	1 or 2 points	●
Barometric Pressure Correction	60.0~112.5kPa/450~850mmHg, manual	●
Salinity Correction	0.0~50.0g/L, manual	●
General Spec.		
Temperature Compensation	0~100°C/32~212°F, manual or automatic	●
Memory	500 data sets, USB communication interface	●
Power Requirements	3× 1.5V AA batteries or DC5V power adapter	●
Dimensions and Weight	W*D*H170× 85 ×30mm, 300g	●
Package Dimension (W*D*H)(mm)	W*D*H440*320*120mm	●
G.W.(kg)	2.2kg	●